STILO 943914

EN Electric massage gun (USB)
PL Pistolet elektryczny do masażu (USB)
CZ Elektrická masážní pistole (USB)
DE Elektrische Massagepistole (USB)
SK Elektrická masážna pištoľ (USB)
FI Sähköinen hierontapistooli (USB)
HR Električni pištolj za masažu (USB)
LT Elektrinis masažo pistoletas (USB)
LV Elektriskā masāžas pistole (USB)

SLO Električna masažna pištola (USB)
ET Elektriline massaažipüstol (USB)
ES Pistola eléctrica para masaje (USB)
UA Електричний масажний пістолет (USB)
BG Електрически пистолет за масаж (USB)
HU Elektromos masszázspisztoly (USB)
MO Pistol electric de masaj (USB)
RO Pistol electric de masaj (USB)

EN USER MANUAL

1. INTRODUCTION

To ensure safety and proper use, read the instructions before using the product for the first time and keep it, as it contains important information.

2. TECHNICAL SPECIFICATIONS

Product name: STILO	
Battery capacity: 2000 mAh	Vibration frequency: 2000 - 3000 rpm
Charging time: approx. 3 hours	Rated power: 25 W
Input voltage: 7,4 V / 14,8 W	Rated voltage: 3,7 V
Weight: 375 g	Dimensions: 14,8 x 10,7 x 5 cm

3. PACKAGE CONTENTS

Before starting use, check whether all the listed elements are included in the packaging.

- gun
- 4 massage pads
- USB cable

It is recommended to clean the product before first use. When throwing away packaging and plastic elements, do it in places intended for this purpose, with environmental protection in mind.

4. PRODUCT DESCRIPTION

A light, handy massage gun, easy to use, suitable for massaging various muscle groups.

4 interchangeable types of massage heads:

flat head for joint/tendon massage

precision head all parts of the body

U - head for vertebral massage

ball head for muscle massage

5. RULES FOR SAFE USE

WARNING!

Read the user manual carefully before use and strictly follow its instructions when using the product. Please keep the user manual as it contains important information.

1. The product is charged using the included charging kit.
2. It is prohibited to use this product while charging.
3. The time of continuous use of the device should not exceed 20 minutes.
4. The product is not waterproof, do not immerse it in water.
5. Do not attempt to unplug the device from the socket or use it with wet hands.
6. To avoid accidents, do not touch the metal part of the plug.
7. Do not destroy, modify, bend or twist the power cord as this may result in damage.
8. Avoid getting the charging port and power cord dirty.
9. Stop charging if the plug or power cord is damaged or becomes hot during charging.
10. Clean the device to avoid dust and dirt accumulation.
11. Do not store the product in the bathroom, humid places, high temperature, as it may cause electric shock or fire.
12. Do not destroy the product, throw it or modify it. If this is damaged, discontinue use.
13. Do not use the device on irritated, allergic or inflamed skin. In this case, contact your doctor.
14. Do not allow children to use the device.
15. Women who are pregnant or during their menstrual cycle should not use the product without first consulting a doctor.
16. Do not use the device under the influence of alcohol or if you feel unwell.
17. Cleaning and maintenance of the device should be performed by adults.
18. Children are not aware of the potential dangers of using this device. Make sure they do not play with the device.
19. This device is not a toy and should be kept out of the reach of children.

6. INSTRUCTIONS FOR USE

1. Charge the device using the original included power cable and a compatible power brick.
2. In case of complete discharge, the full charging cycle takes 3 hours.
3. To turn on the device, press . The device will start in default mode - level 00. Press the power button again to turn off the device.
4. Using the power adjustment buttons, you can adjust the intensity of the massage according to your individual needs
5. To change the head, grab the replaceable tip and gently twist it to pull it out of the hole.
6. When changing the massage head, to improve its operation and avoid seizures, you can use a special oil-like agent, gently rubbing the walls of the hole from which it was removed.
7. To install the head, place it correctly in the hole and then gently push it in until it stays in place.
8. To turn off the device, press  button. The shutdown will be signaled by an acoustic signal.

Caution: The default operating time of the device is 15 minutes. If you want to extend the use, please restart it.

7. MASSAGE TECHNIQUES

Thanks to the use of advanced technology and ergonomic design, the massage gun provides an intense and at the same time relaxing massage that can be adapted to the individual needs of the user. It is perfect for massaging both the upper and lower parts of the body.

Main principles of using a vibrating massager:

1. Massage the selected muscle group according to the pictures (page 29).
2. Choose the appropriate tip for the body part to be massaged.
3. If you encounter a painful area, hold the massager.
4. Remember to breathe properly – deep inhaled and exhaled.

8. REPAIR AND MAINTENANCE:

1. To clean the device, wipe it gently with a dry cloth. It is advisable to use products for the care of electronic devices.
2. In case of damage or failure, please stop using this product and contact SPOKEY service.
3. If the device is not used for a long time, it is recommended to charge it at least 1-2 times a month and store it in a dry and cool place to extend its life.
4. If you do not intend to use the equipment for a long time, store it in the original packaging, out of the reach of children.

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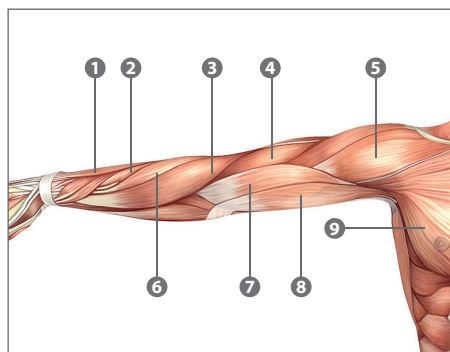
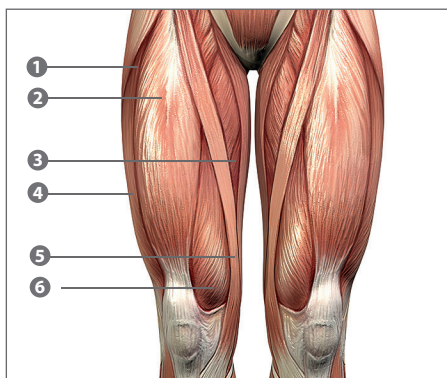
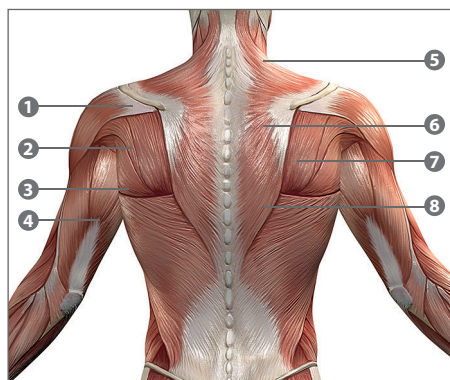
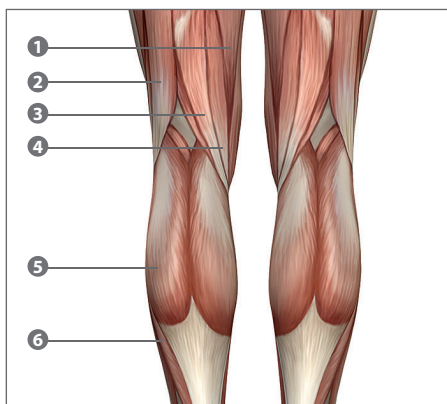
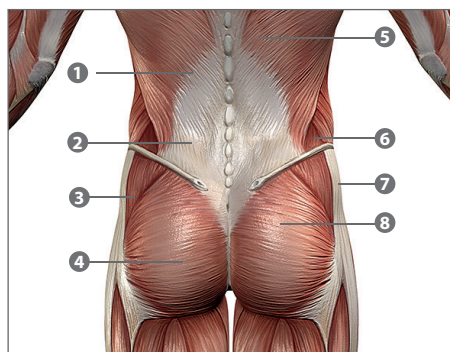
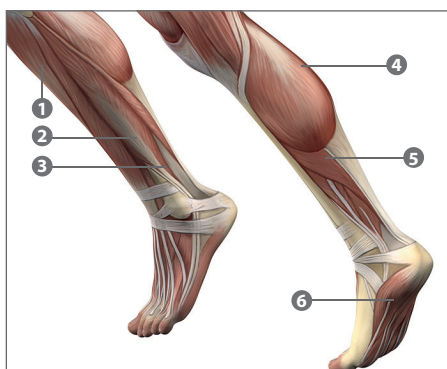
EN DECLARATION OF CONFORMITY WITH EUROPEAN STANDARDS

This device complies with the requirements of the electromagnetic compatibility directive EMC (2014/30/EU), therefore the CE mark has been applied on it and a declaration of conformity with European standards has been issued for it.

This device complies with the requirements of the electromagnetic compatibility directive EMC (2014/30/EU), therefore the CE mark has been applied on it and a declaration of conformity with European standards has been issued for it.

This product meets the requirements of the directive 2011/65/EU and 2015/863 - called the RoHS directive. The purpose of the directive is the approximation of the legislation of the Member States regarding the restriction of use hazardous substances in electrical and electronic equipment and contributing to protect people's health and environmentally friendly waste recycling and disposal of electrical and electronic equipment waste. RoHS directive has been in force since January 3, 2013. New electrical and electronic equipment placed on the market does not contain lead, mercury, cadmium, hexavalent chromium, polybromine diphenyl or polybrominated ether of phenyl ether.

This product meets the requirements of the directive 2011/65/EU and 2015/863 - called the RoHS directive. The purpose of the directive is the approximation of the legislation of the Member States regarding the restriction of use hazardous substances in electrical and electronic equipment and contributing to protect people's health and environmentally friendly waste recycling and disposal of electrical and electronic equipment waste. RoHS directive has been in force since January 3, 2013. New electrical and electronic equipment placed on the market does not contain lead, mercury, cadmium, hexavalent chromium, polybromine diphenyl or polybrominated ether of phenyl ether.

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